


# PressureGuard Turn Select®

Category: Programmable Lateral Rotation Therapy

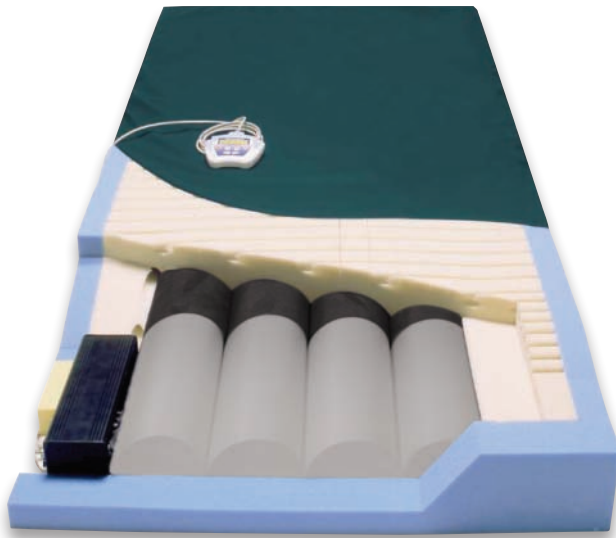

## Customized Lateral Rotation – 24 Hours a Day!

Slowly, gently, and quietly rotates the user laterally up to 30 degrees, with dwell times as short as 1 minute and as long as 2 hours. Hand-held pendant allows convenient programming of the system's four operating modes. Provides caregiver the flexibility to select a standard turning program, or to create a pattern fully customized to the patient's condition.

Ultra-compact, microprocessor controlled pump is integrated directly into the mattress, leaving no bulky pump or external air lines to contend with. Just hour after hour of whisper quiet, trouble-free operation. Maintains proper inflation and support during power loss until power is restored and rotation can be resumed. Includes CPR mode (provides an extra-firm surface for performing CPR with a crash board), program lock-out, and CPR alarm.

### QUICK FACTS

## PressureGuard Turn Select®

#### Primary Use:

- Aggressive Treatment (through Stage IV)

#### Patient Indications:

- Manual turning prohibited by pain
- Can't be positioned off ulcer
- User requires customized turn pattern
- Prevention for those at very high risk (e.g. Braden score of 9 or less)

#### Key Advantages:

- Pendant controlled, no external pump
- Standard or customized roll patterns
- Won't deflate in power outage
- Ultra-quiet
- PressureGuard® design: stability, safety, anti-shearing

Item #	Description
TS7535-29	PressureGuard Turn Select®, 75"L x 35"W x 7"H
* TS8032-29	PressureGuard Turn Select®, 80"L x 32"W x 7"H
TS8035-29	PressureGuard Turn Select®, 80"L x 35"W x 7"H
TS8435-29	PressureGuard Turn Select®, 84"L x 35"W x 7"H

Wider sizes available: See In-Home Styles – page 11.

\* Allow up to 10 working days lead time for this size.

<b>Medicare Approved:</b>	E0277 (Group II)
<b>Weight limit:</b>	350 lbs.
<b>Warranty:</b>	2 years, non-prorated
<b>Mattress Weight:</b>	28 lbs.
<b>UL listed:</b>	Yes

All models are fire barrier version (FBV) as indicated by "-29" suffix, denoting certified conformance to NFPA 101 (Life Safety Code); Cal Tech Bulletin #117; 16 CFR

