



# Geo-Mattress® Series

Though they address different categories of need, the Geo-Mattress Max®, Plus®, Pro, with Wings, Atlas and new *UltraMax* models have much in common. All are protected by covers made from wipe-clean, moisture-proof, anti-microbial fabrics. More importantly, each includes exclusive design innovations that solve problems left unaddressed by other foam mattress replacements:



### Exclusive Geo-Matt® Design

Clinically-proven Geo-Matt® design creates a unique, anti-shearing top surface. 800 individually responsive cells in distinct head, torso, and foot zones enhance therapeutic support and comfort.

(see page 25)


### Hidden Zipper Design

For easier, more thorough cleaning and superior infection control.




### Protective Heel Slope®

Heel Slope® feature redistributes load to pressure-tolerant lower legs and calves. Subtle taper reduces heel pressures while providing complete foot support. Because foot remains level with heart, slope does not facilitate dependent edema.

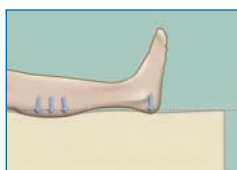

Other surfaces are flat, and can't conform to the shape of the lower leg, especially behind the Achilles tendon.

**Result:** Heel must bear excess load, and skin will break down more easily.

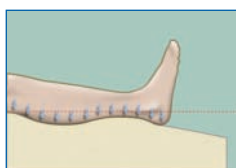

The Geo-Mattress Heel Slope® slightly lowers the heel, allowing the mattress to fill in the space behind the Achilles tendon.

**Result:** Pressure is distributed over the entire lower leg, reducing pressure on the heel.