



Seating

Category: Prevention with Positioning

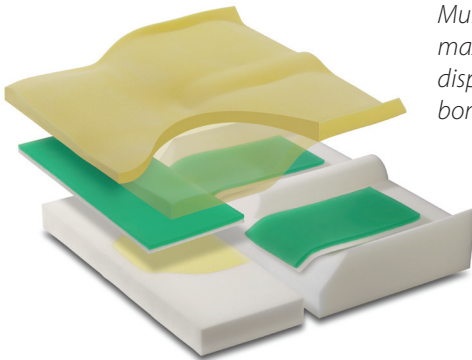

Multi-component construction maximizes absorption and dispersion of pressure from bony prominences

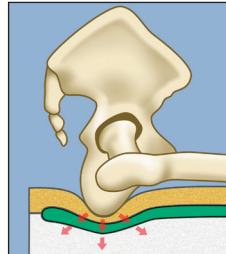

The Equalizer®

The multi-component construction of the Equalizer contoured positioning cushion features an ultra-comfortable viscoelastic “memory foam” top, a contoured foundation layer for proper positioning of the legs, and an extra-soft relief zone for high-risk bony prominences.

Three zones of a proprietary, gel-like soft polymer provide added pressure dispersion under the thighs, greater trochanters, and ischial/tailbone region. Moisture-resistant cover has a rugged nylon bottom with a 2-way stretch fabric top and convenient front carry handle. Built-in sleeve in bottom of cover allows use of molded Sling-Fill™ base (not included- see page 46) for added stability in fabric sling seats; 2-year warranty. Average weight: 4 pounds. Height: 3.5”. Weight limit 250 lbs. (500 lbs. on Bariatric sizes.)

Medicare Code: Skin Protection and Positioning – E2607.
For sizes 22” or wider, use E2608.


EZ-Dish®

A comfortable viscoelastic “memory foam” top molds to the anatomically shaped, round-edged pocket below in the high density bottom. Pocket reduces pressure on ischials and coccyx. Bottom layer is contoured to encourage healthy, upright posture and proper positioning of the legs. Cover consists of anti-slip top surface and non-skid bottom. Built-in sleeve in bottom of cover allows use of molded Sling-Fill base (not included-see page 46); 18-month warranty.

Height: 3”.

Weight limit 250 lbs.

Medicare Code: Skin Protection – E2603.

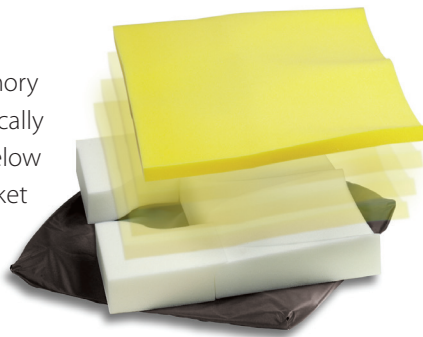

WARNING: Cancer and Reproductive Harm – www.P65Warnings.ca.gov

Item #	Description	Qty/Cs
EQ1616	Equalizer cushion, 16”W X 16”L	1
EQ1616-04	Equalizer cushion, 16”W X 16”L	4
EQ1816	Equalizer cushion, 18”W X 16”L	1
EQ1816-04	Equalizer cushion, 18”W X 16”L	4
EQ1818	Equalizer cushion, 18”W X 18”L	1
EQ1818-04	Equalizer cushion, 18”W X 18”L	4
EQ2018	Equalizer cushion, 20”W X 18”L	1
EQ2018-04	Equalizer cushion, 20”W X 18”L	4

Bariatric sizes:

BEQ2218	Equalizer cushion, 22”W x 18”L	1
BEQ2218-04	Equalizer cushion, 22”W x 18”L	4
BEQ2418	Equalizer cushion, 24”W x 18”L	1
BEQ2418-04	Equalizer cushion, 24”W x 18”L	4

Fire barrier versions can be ordered to conform with 2000 NFPA 101 (Life Safety Code) and Cal TB# 133. To order, add the suffix “-29-1” to the item number for quantity 1, or replace the “-04” suffix with “-29” in the item number for quantity 4.

QUICK FACTS

Use for:

- Pain relief
- Prevention of breakdown
- Assistance with positioning

RISK LEVEL Low Moderate High

Item #	Description	Qty/Cs
EZ1616-96	EZ-Dish w/ anti-slip cover, 16”W x16”L	6
EZ1616-96-1	EZ-Dish w/ anti-slip cover, 16”W x16”L	1
EZ1816-96	EZ-Dish w/ anti-slip cover, 18”W x16”L	6
EZ1816-96-1	EZ-Dish w/ anti-slip cover, 18”W x16”L	1
EZ1818-96	EZ-Dish w/ anti-slip cover, 18”W x18”L	6
EZ1818-96-1	EZ-Dish w/ anti-slip cover, 18”W x18”L	1
EZ2018-96	EZ-Dish w/ anti-slip cover, 20”W x18”L	6
EZ2018-96-1	EZ-Dish w/ anti-slip cover, 20”W x18”L	1


Fire barrier versions can be ordered to conform with 2000 NFPA 101 (Life Safety Code) and Cal TB# 133. To order, replace the “-96” suffix with “-29” in item number.

QUICK FACTS

Use for:

- Pain relief
- Prevention of breakdown
- Assistance with positioning

RISK LEVEL Low Moderate High

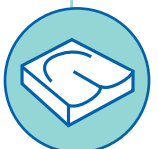

SEATING