

PressureGuard® APM²


Category: Alternating Pressure with Basic Lateral Rotation

Alternating Pressure or Lateral Rotation – at the Flip of a Switch!

With its two distinct therapy modes, the PressureGuard APM² gives you twice the treatment options. Twice the flexibility. Twice the value.

Change instantly and easily from alternating pressure to basic lateral rotation with a single, recessed toggle switch. Requires no disconnecting or exchanging of air lines.

Deluxe control unit includes audible and visual low pressure indicator. This silent, reliable, maintenance-free control unit allows simple adjustment for user comfort. The user experiences subtle but effective alternating pressure or lateral rotation in ten-minute cycles. These dual therapies combine with the patented PressureGuard air/foam design to deliver outstanding tissue load management in any care setting.

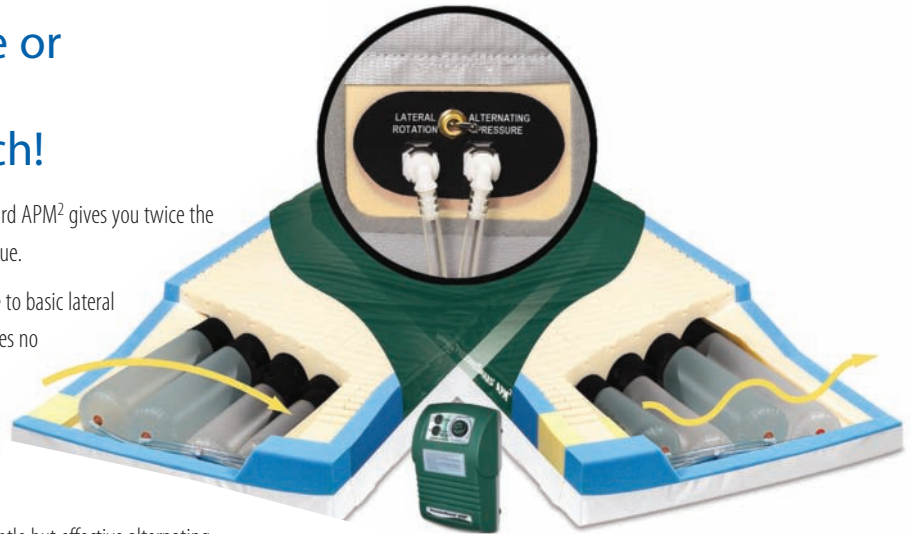

QUICK FACTS

PressureGuard® APM²

Primary Use:

- Aggressive Treatment (through Stage IV)

Patient Indications:

- Can't or won't reposition frequently
- Can be positioned off of ulcer in at least 2 positions
- Prevention for those at high risk for skin breakdown (e.g. Braden score 12 or less)

Key Advantages:

- Two therapies in one mattress
- Single switch operation
- PressureGuard® design: stability, safety, anti-shearing



Deluxe control unit

- Low pressure indicator light
- Selectable audible alarm
- Torque-free ports
- Hospital-grade plug
- Heavy-duty hooks

APM²:

Item #	Description
5875LR-29	PressureGuard® APM ² , w/deluxe control unit, 75"L x 35"W x 7"H
5880LR-29	PressureGuard® APM ² , w/deluxe control unit, 80"L x 35"W x 7"H
5884LR-29	PressureGuard® APM ² , w/deluxe control unit, 84"L x 35"W x 7"H
5800	Deluxe control unit only

Wider sizes available: See In-Home Styles – page 11

APM: Ships with standard control unit, which does not include lateral rotation feature.

Item #	Description
5175-29	PressureGuard® APM, w/standard control unit, 75"L x 35"W x 7"H
5180-29	PressureGuard® APM, w/standard control unit, 80"L x 35"W x 7"H
5184-29	PressureGuard® APM, w/standard control unit, 84"L x 35"W x 7"H
5500	Standard control unit only

Bariatric APM – see page 9.

Medicare Approved: E0277 (Group II)

Weight limit:	350 lbs.
Warranty:	18 month, non pro-rated
Mattress Weight:	20 lbs.
UL listed:	Yes

All models are fire barrier version (FBV) as indicated by "-29" suffix, denoting certified conformance to NFPA 101 (Life Safety Code); Cal Tech Bulletin #117; 16 CFR Parts 1632 and 1633; and BFD IX-11.

