


PressureGuard® APM² Safety Supreme

Category: Multi-Mode Air Therapy with Raised Safety Perimeter




The Ultimate in User Safety, Treatment Flexibility and Skin Integrity.

Combines all the features of the top-selling PressureGuard® APM², plus the added safety benefits of raised perimeter edges. The Safety Supreme incorporates the patented Geo-Mattress with Wings® raised perimeter design to help prompt the user away from the edge of the bed. Provides protection from pressure ulcers, falls, fire, restraints and siderail entrapment – all in one easy-to-use package. **Now with digital control unit.**

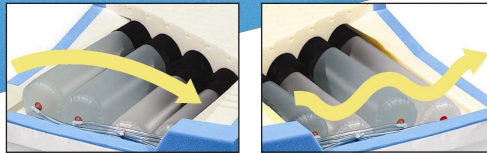

Exclusive Raised Perimeter

incorporates the unique “edge-to-edge” design to help prompt user away from edge of bed without compromising skin integrity.

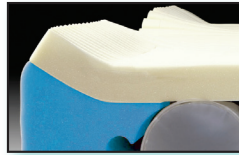

Digital Multi-function Control Unit

Lightweight and quiet, it delivers an unmatched combination of treatment flexibility, ease of use, and dependability.



- “Alternate” mode changes loading in 10-minute cycle; controlled at mattress for choice of lateral rotation or alternating pressure therapy.
- “Float” mode provides powered flotation therapy.
- Timed (20-minute) “Auto Firm” mode for transfers, ADLs, & CPR.
- Simple, indexed comfort adjustment.
- Power Failure indicator light and selectable audible alarm.
- Low Pressure indicator light and selectable audible alarm.

Item #	Description
SAF5875-29	Safety Supreme, 75”L x 35”W x 7”H center, 9”H side
SAF5880-29	Safety Supreme, 80”L x 35”W x 7”H center, 9”H side
SAF5884-29	Safety Supreme, 84”L x 35”W x 7”H center, 9”H side
5900	Digital control unit

Medicare Approved: E0277 (Group II)
Weight limit: 350 lbs.
Warranty: 18 months (all components)
Mattress Weight: Approx. 20 lbs.
UL listed: Yes

Silver3® is a pioneering product developed by medical fabric innovator Dartex Coatings in partnership with recognized wound care leader Smith & Nephew*.

Item #	Description
*84409-29	PressureGuard APM ² Safety Supreme, with Silver3®, 75”L x 35”W x 7”H center, 9”H side
*53667-29	PressureGuard APM ² Safety Supreme, with Silver3®, 80”L x 35”W x 7”H center, 9”H side
*18635-29	PressureGuard APM ² Safety Supreme, with Silver3®, 84”L x 35”W x 7”H center, 9”H side
* Allow up to 10 working days lead time for this version.	

NOW AVAILABLE WITH

Silver3®

proprietary fabric technology

protects surface from stains, odor-causing bacteria, mold and mildew

Silver3® is a registered trademark of Dartex Coatings Ltd • www.dartexcoatings.com • email: sales@dartexcoatings.com • Tel: +44 1159 837 676
 *Uses Smith & Nephew Extruded Films technology • www.snef.co.uk • email: snef@smith-nephew.com • Tel: +44 1430 440 757

