

# Risk Manager®

Category: Bedside Safety Mat

## Minimize the Impact of Falls – On Your Residents and On Your Facility.

Statistics don't lie. Despite the best efforts of the healthcare community, falls by nursing home residents and hospital patients do occur. When they do, the physical, legal, and financial impact can be devastating to all involved.


Span-America's new Risk Manager™ bedside safety mat is designed to help.

The Risk Manager is engineered from an ultra-tough, impact-absorbing lightweight elastomeric gel compound commonly used in expensive orthotics and premium footwear.

Independent laboratory testing (see chart) verifies that the Risk Manager greatly reduces the peak G-forces generated by falls, easily surpassing the impact reduction provided by thicker, bulkier gym mat-style pads. It also outperforms the most commonly purchased specialty fall strip, which relies on the same low-tech foam scrap padding typically found under residential carpeting.

### IMPACT TEST COMPARISON

Drop Height 36"

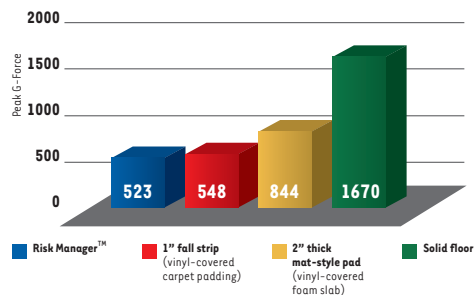

<sup>1</sup>Source: National Technical Systems, Detroit, Michigan. 2008. Results on file.




**Premium elastomeric gel** surface reduces impact from falls. Bonus: Provides welcome comfort for caregivers.



**Nano-silver throughout** protects mat from stains, odor-causing bacteria, mold and mildew. Mat wipes clean with most standard cleaning agents.



**Slide-resistant bottom surface** grips floor, helps keep mat in place. **Textured top** provides added safety.

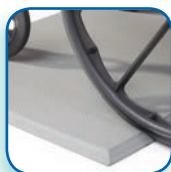

**Low profile and beveled edge** allow easy on and off for wheelchairs and lifts.



**Not a trip hazard**, unlike bulky gym-style mats. Thin, stable surface means staff won't move Risk Manager from beside bed.

Item #	Risk Manager™ Bedside Safety Mat	Qty/Cs
37052	24"W x 70"L x 1"H, standard – single	1/case
94612	24"W x 70"L x 1"H, standard – bulk pack	5/case
8705236	36"W x 70"L x 1"H, XL – single	1/case
9611836	36"W x 70"L x 1"H, XL – bulk pack	5/case

**Warranty:** 1 year  
**Weight:** 6 lbs. (24"), 9.5 lbs. (36")  
**Flammability:** Fire retardant – passes Cal TB 117

**Stain- and oil-resistant**  
**Fluidproof, will not absorb moisture**  
**Non-toxic. Latex-free. Lead-free.**  
**Rolls up for easy transport, storage and shipping.**

### Optional Accessory\*

#### RISK MANAGER™ GRIP STRIPS

For use on highly-polished floors where extra slide resistance may be desired.

\*(Included as standard with XL version.)

**Item # 43398 (3/Package)**  
 (2" wide x 14" long)

