

# PressureGuard CFT®

**Category:** Non-Powered Dynamic Treatment



## Powerful Results — *Without Power*

The revolutionary PressureGuard CFT's patented Constant Force Technology® automatically adjusts a network of interconnected air cylinders and elasticized reservoirs to the appropriate, therapeutic level, regardless of the user's weight or position.

**The CFT is the only advanced, non-powered system of its kind to be recognized for reimbursement as a Group II support surface,** based on head-to-head studies in which it outperformed leading powered systems in ulcer healing.

The CFT is extremely cost effective. Requires virtually no set-up, programming, or maintenance. Wipes clean. Generates no noise or heat and uses no electricity, making it immune to power interruptions.

Use it in-facility to slash rental costs, increase user satisfaction, and improve regulatory compliance on issues from skin care to bed safety. Or, select it as a trouble-free but supremely effective option for at-home users.

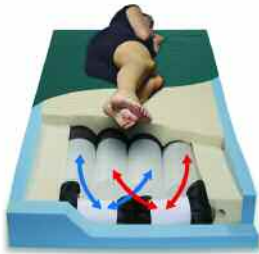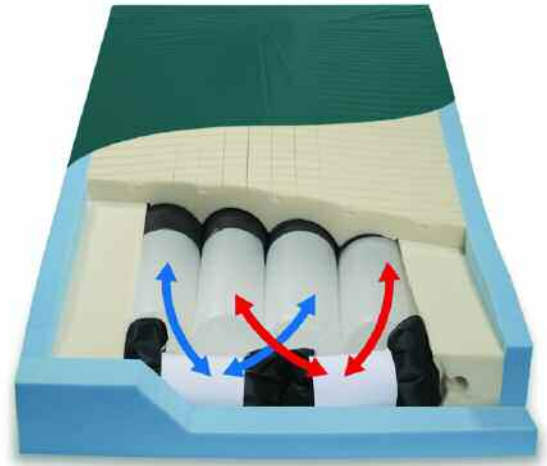

*As user shifts position, reservoirs absorb displaced air to accommodate load and re-equalize pressures.*

Item #	Description	
8075-29	PressureGuard CFT®, 75"L x 35"W x 7"H	
* 8280-29	PressureGuard CFT®, 80"L x 32"W x 7"H	
8080-29	PressureGuard CFT®, 80"L x 35"W x 7"H	
8084-29	PressureGuard CFT®, 84"L x 35"W x 7"H	
* CF8042-29	PressureGuard CFT®, 80"L x 42"W x 7"H	Bariatric
* CF8048-29	PressureGuard CFT®, 80"L x 48"W x 7"H	Bariatric
* CF8054-29	PressureGuard CFT®, 80"L x 54"W x 7"H	Bariatric
* CF8442-29	PressureGuard CFT®, 84"L x 42"W x 7"H	Bariatric
* CF8448-29	PressureGuard CFT®, 84"L x 48"W x 7"H	Bariatric
* CF8454-29	PressureGuard CFT®, 84"L x 54"W x 7"H	Bariatric
SERVICE-00	Air System Renewal Kit (includes hand pump and auto-adjust valve)	

### QUICK FACTS

#### PressureGuard CFT®

##### Primary Use:

- Treatment of any stage pressure ulcer
- Prevention of skin breakdown in the high-risk patient

##### Patient Indications:

- Any stage ulcer
- Can be positioned off of ulcer in at least 2 positions
- Highest risk for skin breakdown, or severely compromised (e.g. Braden score 12 or less)
- History of pressure ulcer(s)
- Repositions infrequently
- Poorly padded bony prominences

##### Key Advantages:

- Proven for both ulcer healing and prevention
- Outperforms powered systems for wound healing\*
- Medicare covered (non-powered Group II)
- No controls to adjust
- No reliance on power source
- PressureGuard® design: stability, safety, anti-shearing

\*Studies available upon request

**Twin/Full/Queen/King sizes available: See In-Home Styles – page 16.**  
**Additional bariatric options available – see page 17.**

\* Allow up to 10 working days lead time for these sizes.

**Medicare Approved: E0373 (Group II)**

**Weight limit: 500 lbs.**

**750 lbs. Bariatric models and Full/Queen/King sizes**

**Warranty: 5 years, non-prorated**

**Mattress Weight: 21 lbs.**

**All models are fire barrier version (FBV)** as indicated by "-29" suffix, denoting certified conformance to NFPA 101 (Life Safety Code); Cal Tech Bulletin #117; 16 CFR Parts 1632 and 1633; and BFD IX-11.